

Rejuvenate your confidence as

JUVESKIN

CROSS-LINKED HA DERMAL FILLER



JUVESKIN

CROSS-LINKED HA DERMAL FILLER



JUVESKIN HA Dermal Filler with lidocaine represents a new trend in aesthetic enhancement. JUVESKIN HA Dermal Filler produced through a safe & stable process using high purity and high molecular weight hyaluronic acid have already been recognized for their excellence in the market of many countries. You can get the most natural result by using the optimal version for each area and condition

3C Technology of JUVESKIN

Conservative Cross-linking system

The 14-day process, low-temperature crosslinking process slowly mixes the crosslinked molecular weight and free HA to ensure uniform texture with 1-2 micron particles and preserves viscoelastic properties that allow for better tissue integration.

Certified safe raw material

We use high-purity HA whose safety has been fully verified worldwide.

Conservative Cross-linking system

Minimize differences between production batches through strict quality control

Injection Force Test

JUVESKIN has a gel-like consistency with a regular, dense, and firm lattice structure as monophasic structure.

This allows for stability, consistency and gentle injection protocols.




Normal force (N)	SOFT	MEDIUM	HARD	LIPS
	24	14	8.3	21





Natural filling with
soft viscoelasticity

Specification

JUVESKIN SOFT	1ml
HA concentration	24mg
Lidocaine HCl	3mg/ml
G'	85
Recommended needle	30G
Complex Viscosity(η^*)	142
Injection Depth	Superficial Dermis
Normal force (N)	24

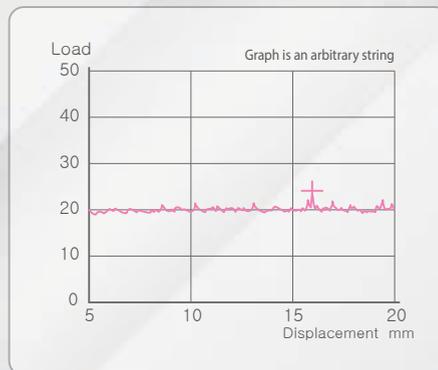

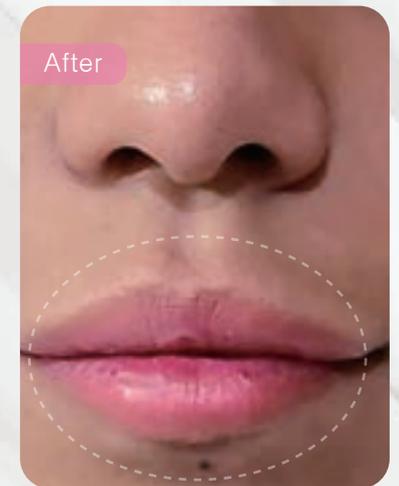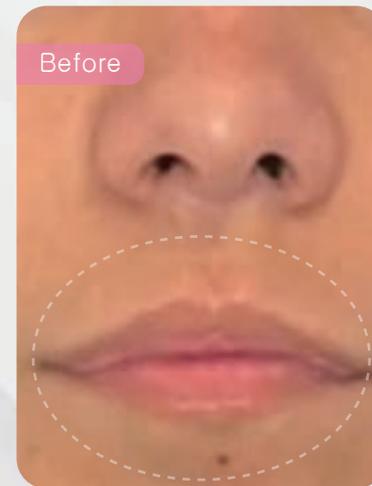








Natural filling
with firm elasticity

Specification

JUVESKIN MEDIUM	1ml
HA concentration	24mg
Lidocaine HCl	3mg/ml
G'	285
Recommended needle	27G
Complex Viscosity(η^*)	460
Injection Depth	Mid/Deep Dermis, Subcutaneous
Normal force (N)	14

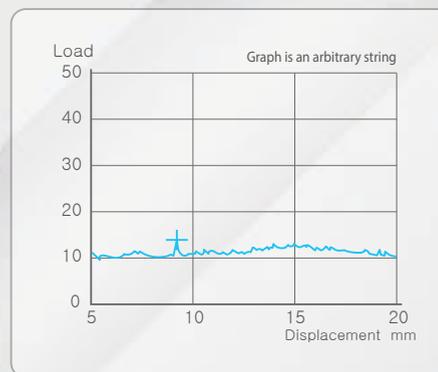

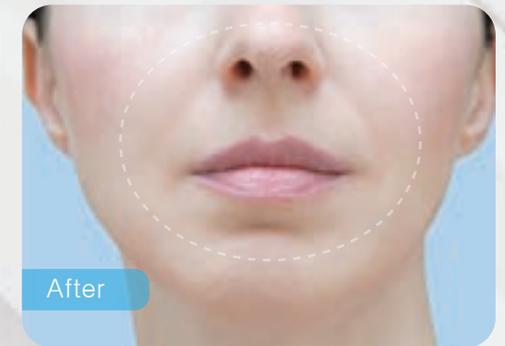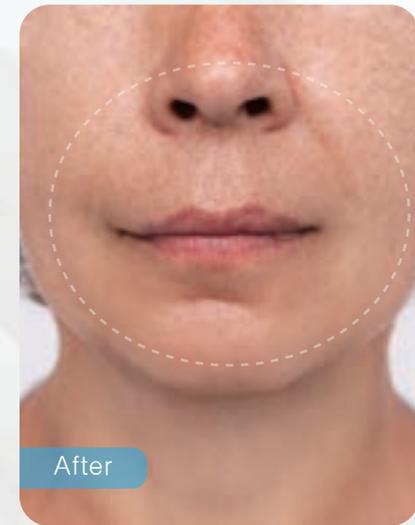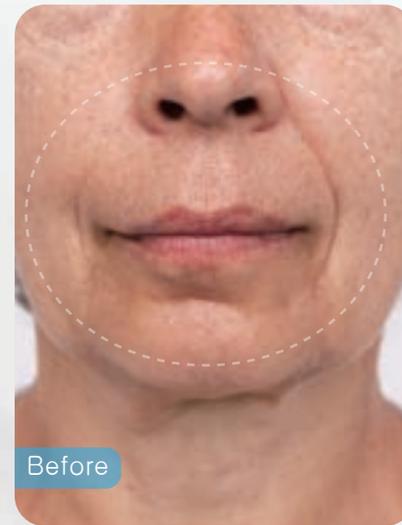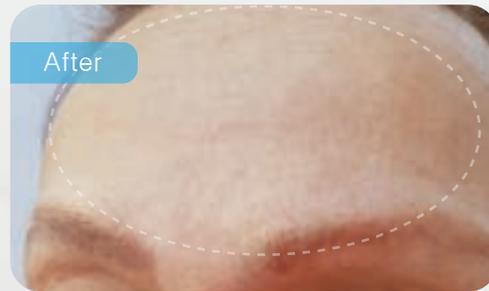






Natural filling of
desired shape with
solid viscoelasticity

Specification

JUVESKIN HARD	1ml
HA concentration	24mg
Lidocaine HCl	3mg/ml
G'	470
Recommended needle	25G
Complex Viscosity(η^*)	762
Injection Depth	Deep Subcutaneous
Normal force (N)	8.3

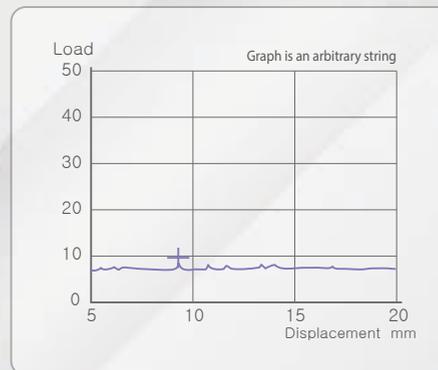

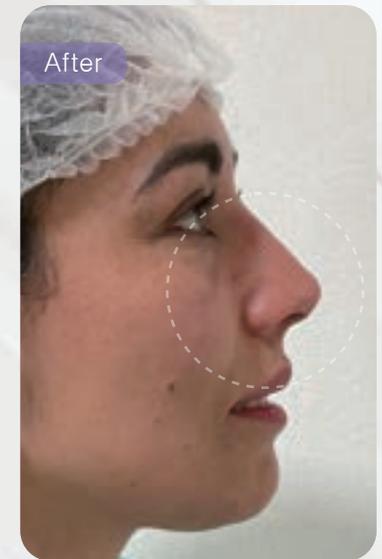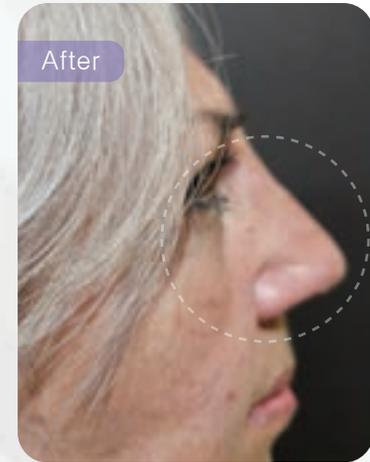



Natural coloring
with vitamin B12

Specification

JUVESKIN LIPS	1ml
HA concentration	24mg
Lidocaine HCl	3mg/ml
G'	85
Recommended needle	30G
Complex Viscosity(η^*)	170
Injection Depth	Superficial Dermis
Normal force (N)	21





JUVESKIN

HA VOLUME AUGMENTATION



The highly purified viscoelastic body filler JUVESKIN BODY HA Filler is the best solution for replenishing body volume and professional body contouring.

JUVESKIN BODY HA FILLER **restores volume to the body** in areas like breasts, lower legs, and buttocks. It is also efficient in combating cellulite and uneven skin surface due to liposuction procedures.

3C Technology of JUVESKIN

Conservative Cross-linking system

The 14-day process, low-temperature crosslinking process slowly mixes the crosslinked molecular weight and free HA to ensure uniform texture with 1–2 micron particles and preserves viscoelastic properties that allow for better tissue integration.

Certified safe raw material

We use high-purity HA whose safety has been fully verified worldwide.

Conservative Cross-linking system

Minimize differences between production batches through strict quality control

Hyaluronic acid in JUVESKIN Body HA Filler has the highest level of purification for maximum safety, leaving minimal residue of cross-linking agent, endotoxins, and proteins.






Benefits

- JUVESKIN BODY HA Filler cross-linking that ensures a longer effect
- Contains highly-purified stabilized hyaluronic acid
- Has higher viscosity that prevents deformation
- Smooths the skin and combats cellulite
- Suitable for patients, who are not intended to go through surgery, but desire to enlarge volume in areas like breasts and buttocks
- Hyaluronic acid gel does not contain animal components, eliminating the risk of allergic reactions
- General anesthesia is not required
- No incisions or scars
- Quick recovery period (compared with plastic surgery)

Indications

- Asymmetrical breasts and buttocks
- Uneven skin surface, including cellulite
- Volume loss in breasts and buttocks
- Contour correction of genitals (enlargement of labia majora, penis)
- Subcutaneous tissue loss due to liposuction
- Funnel chest
- Asymmetric calves
- Soft tissue defects related to surgery and trauma

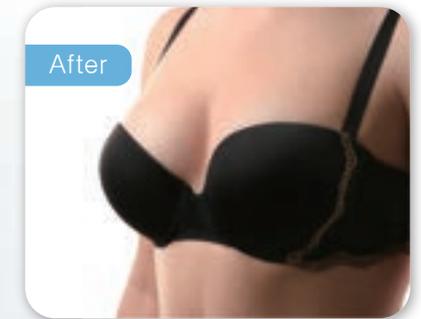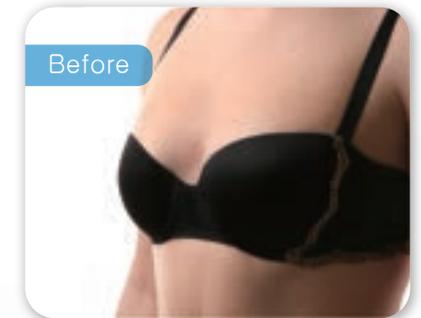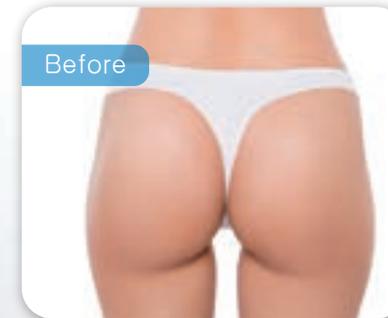


NOTICE

The majority of our products are designed for professional use by doctors and specialists, by purchasing said product(s) you agree to supply on request formal proof of medical/professional training. Said product(s) are not intended for personal or home use.

Please note this product has contraindications for use. It is recommended that all Healthcare professionals familiarise themselves with these contraindications and precautions by reading the instructions provided.

Unique Technology

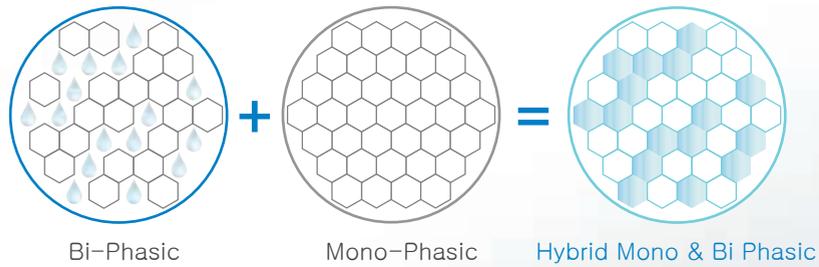

- High elasticity
- Smooth injection
- High quality and safety
- Natural volume

Body Contour Filler "JUVESKIN BODYFILL 80"

Characteristics

The experimental group showed better physical properties than the control group.

- High cross-linked HA gel content compared to similar products
- High viscoelasticity of cross-linked HA gel suitable for body filler
- Soft and low injection force compared to high viscoelasticity

Products	HA Content (mg/ml)	G' (Pa) Storage modulus	η^* (Pa.S) Complex viscosity	G''(Pa) Loss modulus	Tan Delta (G''/G')	Injection Force (N)
JUVESKIN	24	470.6	430.0	92.8	0.20	9.14
Control A	20	395.4	505.5	43.7	0.11	12.42
Control B	20	212.4	435.3	41.1	0.19	11.34
Control C	20	163.8	205.6	28.2	0.17	10.34





Rejuvenate your confidence as

JUVESKIN

CROSS-LINKED HA DERMAL FILLER



DFKBIOLAB INC.

ADDRESS. #522, 272 Sunhwagung-ro, Namyangju-si, Gyeonggi-do, 12106, Republic of Korea

TEL. +82-2-6225-7988 (7985,7984) **E-mail.** dfkbiolab@sbodyline.com **HOME PAGE.** www.sbodyline.com